

- [Recording Instructions and Requirements](#)
 - [Pre-filming requirements:](#)
 - [During-filming requirements:](#)

Recording Instructions and Requirements

For the recording to be properly recognized by the AI model. A few things need to be marked correctly such that it sees the desired points correctly for appropriate post-processing.


The model is designed to detect: 4 dots on the skin that represent the 4 corners of a square with sides of 1 or 2 cm. A filament with 3 separate color zones, in order to distinguish 6 points for bending.

Pre-filming requirements:

- Static clean white lighting for the subject.
- Stationary camera.
- Paint the dots a bright opaque green. e.g. [Posca](#) has some paint-type markers, it was what was used for the training. e.g. Bright green dots on skin, in a white-lit setting:



- Paint the filament a bright opaque [colour] for 3 defined segments. e.g. Example of a filament painted white and dark blue:

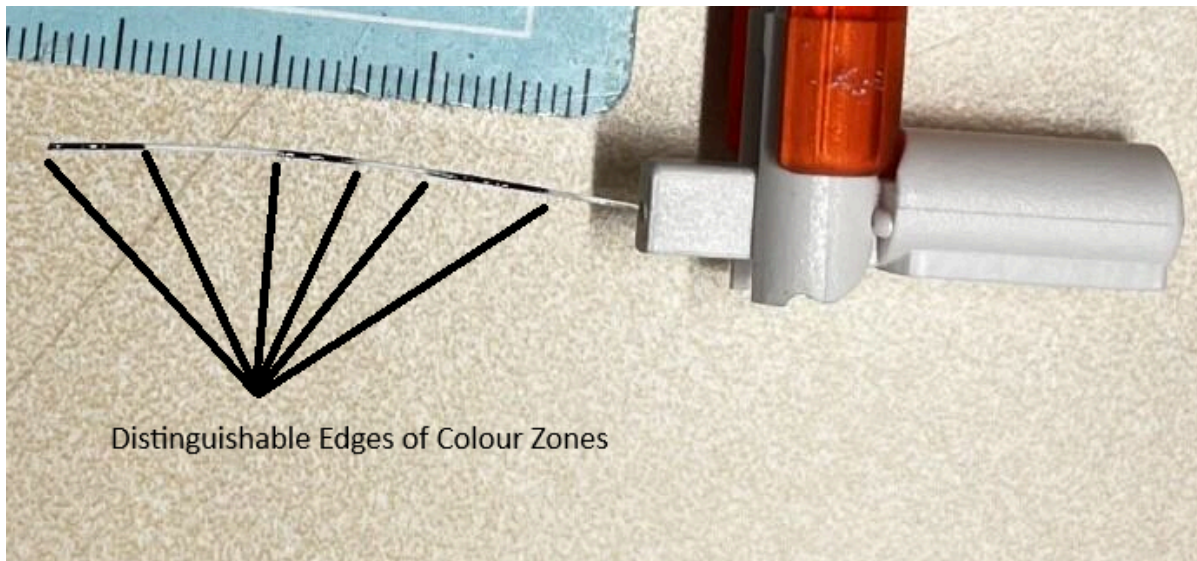

- Clean any old dots/blemishes on the skin.

During-filming requirements:

- For synchronization with neuron data; Start the video with 5 touches, roughly 1 per second, at the hotspot.
- No filaments elsewhere within frame during recording.
- Double check recording for bending visibility Image 1: Region of interest, the filament, is very blurry. Image 2: Bend is perpendicular to the cameras reference. Image 3: Good image, not blurry, very visible bend.

